

Combination Blanket/Fluid Warmer

P-2147 P-2150





MN-39930

REV.04 9/23

ΕN

Manufacturer's Information

Copyright	© Copyright 9/23 by Pedigo Products, Inc.		
	All rights reserved.		
	This manual or any portion thereof may not be reproduced or used in any manner whatsoever without the express written permission of Pedigo Products, Inc.		
Trademarks	All trademarks referenced in this documentation are the property of their respective owners.		
Manufacturer	Pedigo Products, Inc.		
	4000 S.E. Columbia Way		
	Vancouver, WA 98661, USA		
	1-800-822-3501		
	www.pedigo-usa.com		
Original instructions	The content in this manual is written in American English.		



Thank you for your Purchase!

This warmer has been thoroughly tested and inspected to ensure only the highest quality is provided. We supply the most durable, convenient, efficient and safe warming equipment on the market. All warmers are manufactured and fully inspected in the USA with a commitment to quality.

Pedigo Repair Service

Call

Call 262-509-6517 to reach our service call center for immediate access to local authorized service agencies. This service is provided exclusively for Pedigo Product, Inc. equipment and is available throughout the United States.



This page is intentionally left blank.



Foreword

Manufacturer's Information
Foreword3Thank you for your Purchase!3Pedigo Repair Service.3
Safety7The Meaning of Signal Words7Safety Precautions.8
Labels11Label Locations
Specifications13Specification Information
Installation21How to Unpack the Warmer.1How to Install the Warmer.23
Operation27Preparing the Warmer for First Use27How to Turn On and Turn Off the Warmer28How to Operate the Warmer30How to Change the Temperature Scale32How to Change the Sound Settings32How to Lock and Unlock the Controller33How to Download Temperature Data from the Warmer34
Maintenance37Maintenance Schedule
Troubleshooting43What to do if a Power Interruption Occurs43What to do if the Alarm Indicator Light Flashes.44Error Codes45
Electrical Information 49
Warranty 53 Warranty



This page is intentionally left blank.



The Meaning of Signal Words

This manual contains signal words where needed. These signal words must be obeyed to reduce the risk of death, personal injury, or equipment damage. The meaning of these signal words is explained below.



DANGER

Danger indicates a hazardous situation which, if not avoided, will result in serious injury or death.



WARNING

Warning indicates a hazardous situation which, if not avoided, could result in serious injury or death.



CAUTION

Caution indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.



Notice indicates a situation which, if not avoided, could result in property damage.



NOTE: Note indicates additional information that is important to a concept or procedure.



Safety Precautions

Before you begin	Read and understand all instructions in this manual.				
Electrical precautions	Obey these electrical precautions when using the warmer:				
	Connect the warmer to a properly grounded outlet. Do not use the warmer if it is not properly grounded. Consult an electrician if there is any doubt that the outlet used is properly grounded.				
	Keep the cord away from hot surfaces.				
	Do not attempt to service the warmer or its cord and plug, when plugged in.				
	Do not operate the warmer if it has a damaged cord or plug.				
	Do not immerse the cord or plug in water.				
	Do not let the cord hang over the edge of a table or counter.				
	Do not use an extension cord.				
Usage precautions	Obey these usage precautions when using the warmer:				
	 Only use this warmer for its intended use of warming medical solution bags, bottles, and/or blankets. Follow facility and solution manufacturer guidance regarding warming temperature for any item to be placed in the warmer. 				
	Do not use this warmer for warming blood or blood products.				
	Do not cover or block any of the openings of this warmer.				
	Do not use this warmer in a wet location.				
	Only clean the warmer when the power cord is unplugged.				
	Do not use corrosive chemicals when cleaning the warmer.				
	Do not use the warmer cavity for storage.				
	Do not remove exterior panels from the warmer or attempt repairs. The warmer has no user-serviceable internal components. Only perform routine cleaning and maintenance procedures specifically described in this manual. Inspection and servicing of internal components must only be performed by qualified service personnel.				
	 Only a qualified Pedigo service representative may make modifications to the warmer. Modifications to the warmer could be hazardous to users and patients. 				
Operator training	All personnel using the warmer must have proper operator training. Before using the warmer:				
	Read and understand the operating instructions contained in all the documentation delivered with the warmer.				
	Know the location and proper use of all controls.				
	Keep this manual and all supplied instructions, diagrams, schematics, parts lists, notices, and labels with the warmer if the warmer is sold or moved to another location.				





Operator qualifications	Only trained personnel with the following operator qualifications are permitted to use the warmer:				
	Have received proper instruction on how to use the warmer.				
	Are familiar with the purpose, limitations, and associated hazards of the warmer.				
	The warmer must not be used by:				
	People impaired by drugs or alcohol.				
Condition of warmer	Only use the warmer when:				
	 All controls operate correctly. 				
	The warmer is installed correctly.				
	The warmer is clean.				
	The warmer labels are legible.				
Servicing the warmer	 Obey precautions in the manual, on tags, and on labels attached to or shipped with the warmer. 				
	 Only trained personnel are permitted to service or repair the warmer. Repairs that are not performed by a trained technician, or the use of non-factory parts, will void the warranty and relieve all liability. 				
	 Any troubleshooting guides and components views included with this manual are for reference only and are intended for use by qualified and trained service technicians. 				
	To prevent serious injury, death or property damage, have the warmer inspected and serviced at least every twelve (12) months by a trained technician.				
	Contact Pedigo for the authorized service partner in your area.				
Incident notice	Any serious incident that has occurred in relation to the warmer should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.				





This page is intentionally left blank.



Label Locations



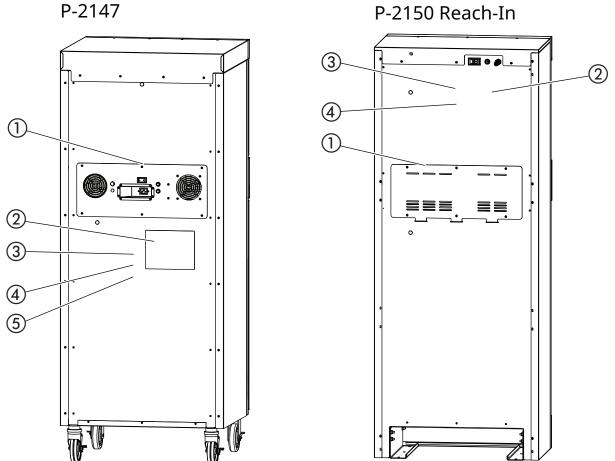

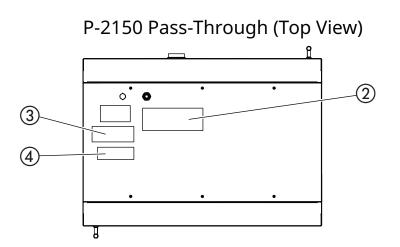

ABEL



EC-PHD-009217

1	Social So	Security seal
2	NAMENION ACCESS BIOLINE DE MARGE BY CALILIPIE SERVER TROVINCIAS ONT. ACCESS DOI TIME MESSARE SERVER SCULSIVIALMAN ACCESS DOI TIME MESSARE STRATEGING CALILIPIE SERVER TROVINCIAS ONT. ACCESS DOI TIME MESSARE SCULSIVIALMAN ACCESS DOI TIME MESSARE STRATEGING CALILIPIE ACCESS DOI TIME ACCESS ON THE ACCESS ON THE COLUMN ACCESS ON THE ACCESS ON THE COLUMN ACCESS ON THE ACCESS ON THE COLUMN ACCESS ON THE ACCESS ON THE ACCESS DOI TIME ACCESS ON THE ACCE	Electrical warning label
3	Distributed by: Adde in Pecipication Adde in Www.pedigo-usa.com Pedigo Products, Inc. 4000 SE Columbia Way Vancouver, WA 98661 800-822-3501 For Service or Emergencies Call 262-509-6517	Pedigo contact information
4	QUALITY INSPECTED	Quality inspected
5	WARNING: This product can expose you to chemicals including Formaldehyde, which are known to the state of California to cause cancer. For more information, go to www.p65warnings.ca.gov. AVERTISSEMENT: Ce produit peut vous exposer à des substances, helles que la formaldéhyde, déclarées cancérigines par l'Etat de Californie. Pour plus d'information, aller à www.p65warnings.ca.gov. ADVERTENCIA: Este product to peut exposer a substancia squincias como el formaldéhido que, según el conocimiento del estado de California, provoca cáncer. Para obtener más información, visite www.p65warning.ca.gov LA47601	WARNING This product can expose you to chemicals including Formaldehyde, which are known to the state of California to cause cancer. For more information, go to www.p65warning.ca.gov.
_		WARNING Hot surface (label installed on the interior cavity wall)



Specification Information


	Model	P-2147	P-2150		
	Intended use Bottom cavity—warming blankets or irrigation fluids Top/TriMode cavity—warming blankets, irrigation or injection flui				
		Net: 457 lb (207 kg)	Net: 336 lb (152 kg)		
	Weight	Ship: 528 lb (239 kg)	Ship: 407 lb (185 kg)		
	Storage cavity capacity	Top/TriMode cavity: 24 1-liter bottles or 20 1-liter bags, or 9–12 blankets (3.9 ft ³)	Top/TriMode cavity: 30 1-liter bottles or 24 1-liter bags, or 10–14 blankets (5.4 ft ³)		
		Bottom cavity: 24–30 blankets (9.4 ft ³) or 30 bottles	Bottom cavity: 30–40 blankets (13 ft ³) or 40 bottles		
			fluids (IRR): (32°C to 66°C)		
	Temperature range	Injection fluids (INJ): 90°F to 104°F (32°C to 40°C)			
		Blanket: 90°F to 160°F (32°C to 71°C)			
	Clearance requirements	4" (102mm) from rear 1" (25mm) from top 2" (51mm) from sides 1/2" (13mm) from bottom	2" (102mm) from rear 1" (25mm) from top 1" (51mm) from sides NOTE: See topic <i>Dimension Drawings</i> for P-2150 Built-In requirements.		
Transportation and	Ambient temperature range of -40°F to +159°F (-40°C to +70°C).				
storage conditions	Relative humidity range of 10% to 95%, non-condensing.				
	Atmospheric pressure range of 7.25 psi to 15.37 psi (50kPa to 106kPa).				
Operating conditions	The warmer must acclimate to the room temperature it will be placed in—24 hours is recommended.				
	 The recommended environmental temperature range is 60°F to 90°F (15°C to 32°C). 				
	The recommended relative humidity is above 20%, non-condensing.				
Standards for electrical equipment	 Laboratory equipment listed by Underwriters Laboratories with respect to electrical shock, fire, and mechanical hazards only, in accordance with UL 61010-1 and CAN/CSA C22.2 No. 61010-1. 				
	 Grounding reliability can only be achieved when the appliance is connected to an equivalent receptacle marked "Hospital Grade." 				
	 EN61326-1; EN61000-3-2; EN61000-3-3 				
	 IP-XO ordinary 				
		,			




PEDIGO

Dimension Drawings

P-2147 shown with feet

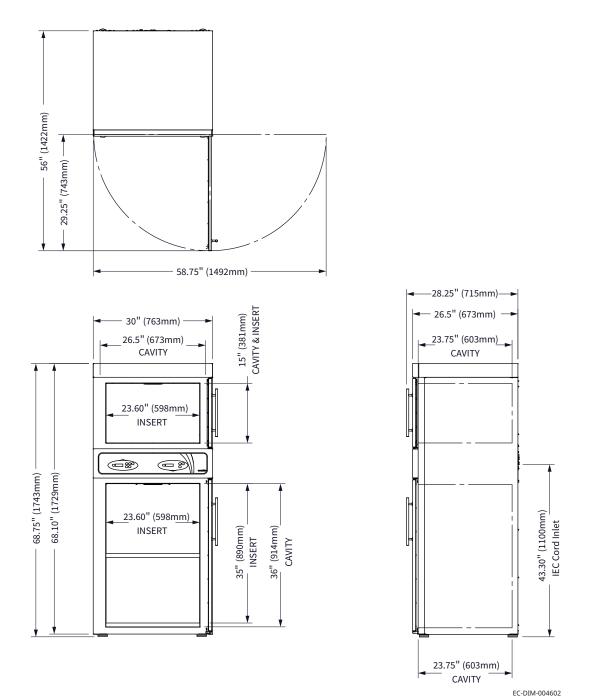



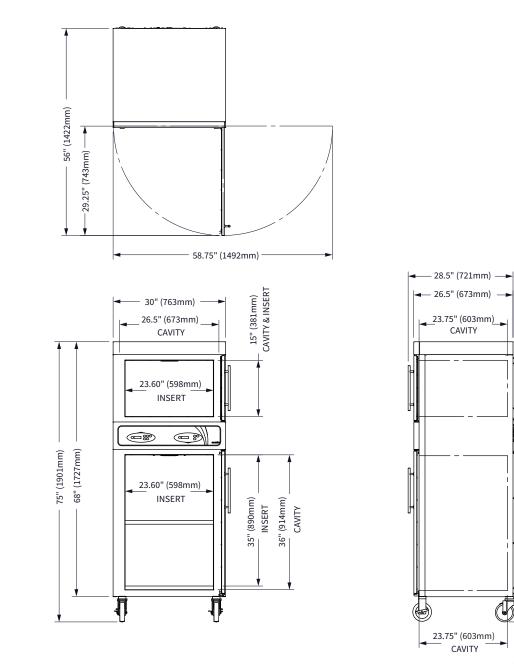


P-2147 shown with casters

EC-DIM-004605

49.5" (1256mm) IEC Cord Inlet



P-2150 Reach-In

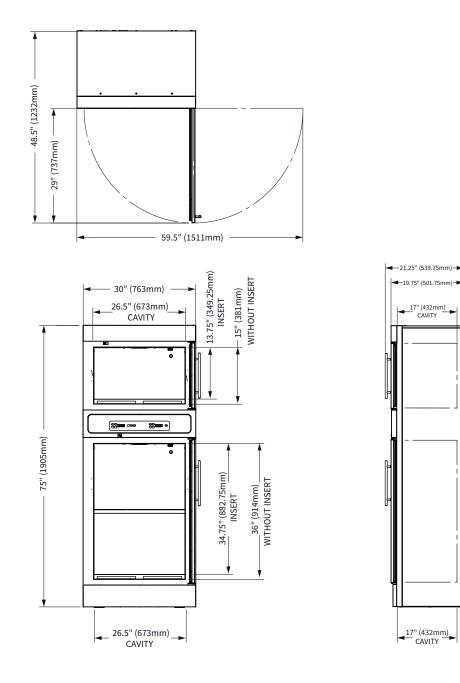




72.75" (1847.75mm) IEC Cord Inlet


EC-DIM-008782


P-2150 Pass-Through

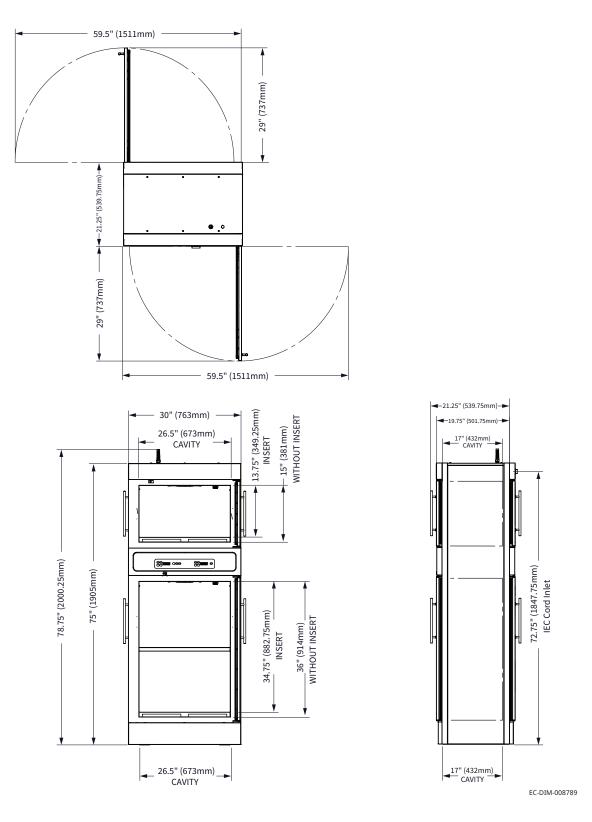



P-2150 Built-In

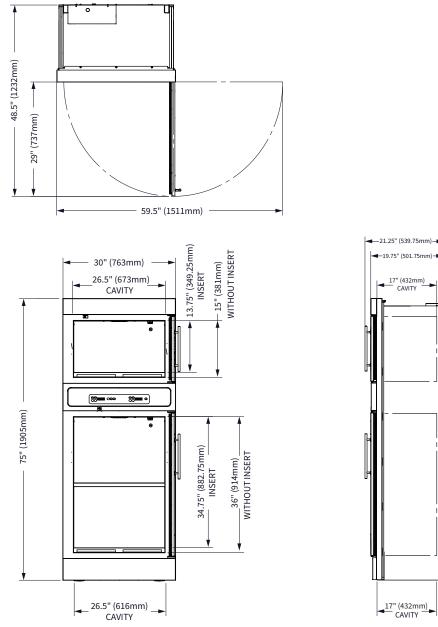

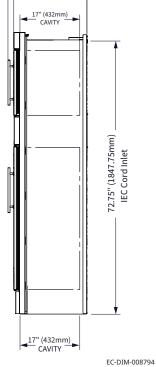



P-2150 Built-In Cutout Specifications and Inlet Locations



CAUTION: Fire hazard. Supply wires must be rated for a minimum wire temperature of 80°C.

NOTICE Built-in units must be hardwired to the facility mains with a 20A breaker switch.

Breaker*	Wire Size**	Wire Temperature
20A	12 AWG (1mm ²)	80°C (minimum)

*Electrical connections must meet applicable codes. **IEC 60335 minimum.


Cutout Specifications			
Height	74.5" (1892.25mm)		
Width	29.25" (743mm)		
Depth	17" (432mm)		
Top of unit clearance	12" (304.75mm)		

Inlet Options				
1	Back of warmer			
2	Top of warmer			
3	Left side of warmer			
_	—			

Continued on next page



Continued from previous page

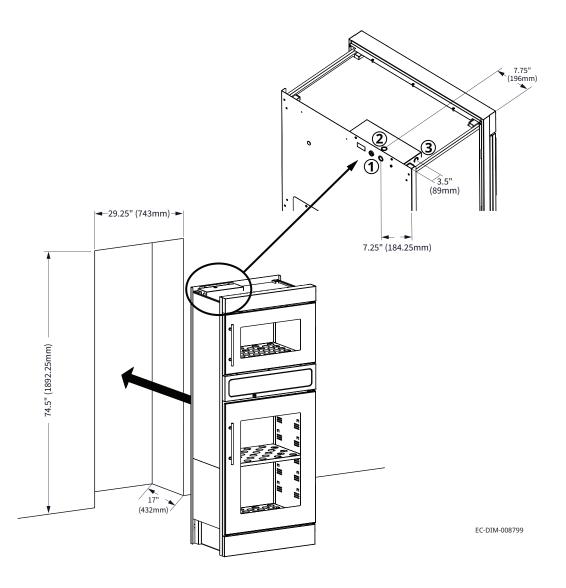



How to Unpack the Warmer



WARNING: Crushing hazard.



An unstable warmer can lead to a crushing hazard. Observe your

facility's best practices for moving large equipment.

Before you begin

Make sure you have:

- Reviewed the warmer specifications.
- An appropriate lifting device and enough personnel to safely move and position the weight of the warmer.
- Cutting tools to remove the packaging.


Unpack the warmer

To unpack the warmer, do the following.

Step	Action			
1.	Remove the box or crate. Save all packing materials for inspection by the carrier.			
	NOTE: Examine the warmer for damage. If the warmer has been damaged, do not use the warmer until it has been inspected by an authorized service provider. Contact your carrier or customer service.			
2.	Cut and remove the retaining straps and plastic wrap.			
3.	Remove the warmer from the pallet using an appropriate lifting device.			
4.	Remove the paperwork from the cavity.			

Continued on next page



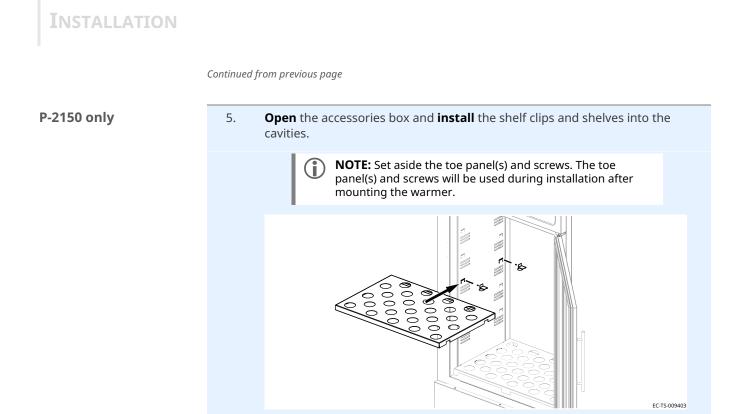

Result

The warmer is now unpacked.



How to Install the Warmer



WARNING: Tipping hazard. Mount the P-2150 to the wall or floor before use.

Before you begin	Make sure you have:				
	An appropriate lifting device and enough personnel to safely move and position the weight of the warmer.				
Requirements	The warmer must be installed on a level surface.				
	The warmer must not be installed in any area where it may be affected by dripping water, high temperature, or any other severely adverse conditions.				
Voltages					

Model	v	Ph	Hz	kW
P-2147	120	1	50/60	1.3
P-2150	120	1	50/60	1.2

Position the warmer

To position the warmer, do the following.

Step	Action
1.	Make sure that:
	 The location where the warmer is being installed is rated to support the weight of the warmer;
	The warmer is within five feet of the appropriate electrical outlet;
	 You follow the warmer clearance guidelines shown below. Refer to topic Specification Information for clearance requirements.
	 You reference the built-in specification page for cutout dimensions (P-2150 only).
2.	Move the warmer to the installation location and onto the final resting surface.
3.	Lock the casters, if equipped.

The warmer is now correctly positioned.



Mounting the warmer to wall or floor (P-2150 only)


To mount the warmer to the wall or floor, do the following.




WARNING: Tipping hazard. Mount the P-2150 to the wall or floor before use.

NOTE: Reference OPM-0562-19 for California seismic conformance.

4. If wall mounting, **bolt** the bracket to the wall using facility approved fasteners (not included).

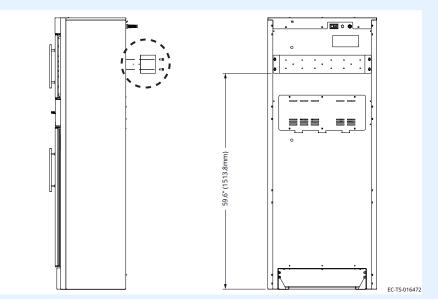


If floor mounting, **bolt** the warmer to the floor using facility approved fasteners (not included) or contact Customer Service for optional floor mounting hardware.

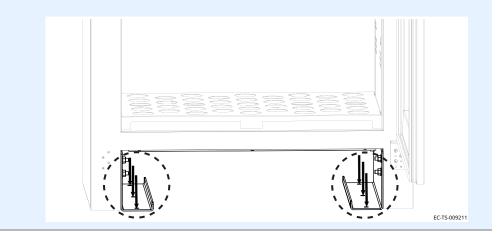

Continued on next page



Continued from previous page

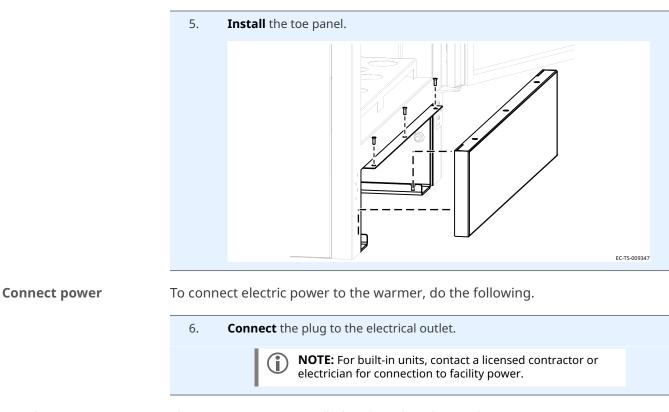

Result

The warmer is now installed and ready to be used.



This page is intentionally left blank.



Preparing the Warmer for First Use

Before you begin	L	CAUTION: Burn hazard. Allow the warmer to cool before cleaning.				
	ľ	NOTICE Do not use: abrasive cleaning compounds. chloride based cleaners. commercial or household cleaners containing ammonia.				
Approved cleaning agents		ი greater isopropyl alcohol Sani-Cloth® HB				
Procedure	To prep	pare the warmer for first use, do the following.				
	Step Action					
	1.	Make sure that the warmer is turned off and cool.				
	2.	Wipe the outside of the warmer:				
		with a stainless steel cleaner, if stainless steel panels.				
		with an approved cleaning agent, if painted panels.				
	3.	Clean the interior of the warmer with a damp cloth or approved cleaning agent.				
	4.	Dry the interior and exterior of the warmer with a clean, lint-free cloth.				
		Leave the door open until the interior of the warmer completely dries.				
	5.	Clean each side of the window pane with glass cleaner or distilled vinegar.				

Result

The warmer is now ready for use.




How to Turn On and Turn Off the Warmer

Before you begin

The warmer must be connected to electric power.

Turning on the warmer

To turn on the warmer, do the following.

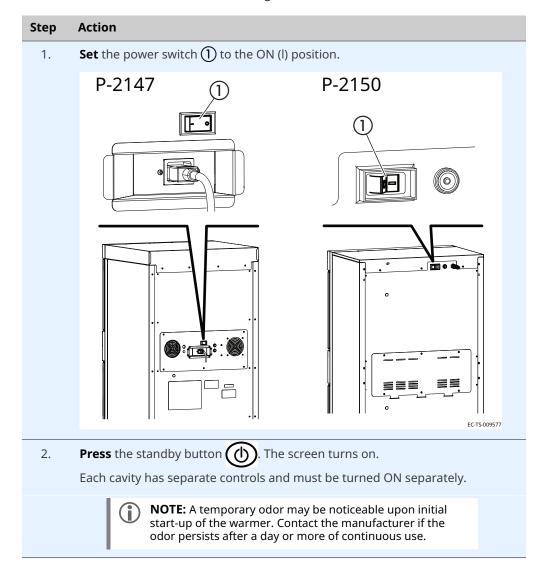

The warmer is now on.

Continued on next page



Continued from previous page

Turning off the warmer

To turn off the warmer, do the following.

Press and hold the standby button until the screen turns off.
 Each cavity has separate controls and must be turned OFF separately.

The warmer is now off.




How to Operate the Warmer

	WARNING: Personal injury hazard.				
	Do not operate the warmer in the presence of a flammable anesthetic mixture with air, oxygen, or nitrous oxide; in oxygen-enriched environments; or in any other potentially explosive environment.				
	WARNING: Personal injury hazard.				
	Make sure all parts shipped with the warmer (trays, shelves, baskets, false bottoms, etc.) are installed to prevent heat degradation of medical products placed within.				
	Do not operate the warmer at a temperature that exceeds the instructions for use of the medical product warmed within.				
	To reduce the risk of accidental set-point adjustment, injectable medical fluid products should only be warmed in a cavity that restricts the maximum set-point.				
	Prior to use of a medical fluid product, verify the temperature is in accordance with the fluid manufacturers instructions for use. Do not use any fluids that are above the recommended temperature.				
NOTI	CE Do not overload the cavity.				
	Allow 1" (25mm) between medical products to be warmed and the interior walls, see topic <i>Specification Information</i> for the storage capacity.				

Procedure

To operate the warmer, do the following.

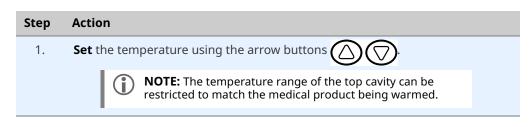

Continued on next page



Continued from previous page

	Step	Action
	1.	Select the operation mode.
		Press the irrigation button to restrict the maximum set-point to 150°F (66°C).
		• Press the injection button $vert$ to restrict the maximum set-point to 104°F (40°C).
		If irrigation or injection are selected, the operation mode will illuminate 🔗.
		 Press the blanket button Sector to restrict the maximum set-point to 160°F (71°C).
		The blanket operation mode icon will not illuminate.
		NOTE: The maximum temperature in the lower cavity is 160°F (71°C). It is recommended that the temperature is locked when warming irrigation fluids in the lower cavity to prevent accidental changes.
	2.	Open the door and load the fluid or blankets into the warmer.
		Close the door.
		NOTICE Do not warm items containing plastic, rubber, metal snaps, studs, hooks, etc.
During the warming process	3.	Press the temperature recall button to view the measured cavity temperature.
		The measured cavity temperature will display for five seconds. Then, the set-point temperature will display.
Result	The flui	ds or blankets are now warming.



How to Change the Temperature Scale

Before you begin	The warmer must be connected to electric power and the screen turned					
Procedure	To change the temperature scale from °F to °C and vice versa, do the following.					
	Step Action					
	1.	Press and hold the standby button with the screen turns off, then release the button.				
	2.	Press and hold the temperature recall button Until the temperature digits disappears and only the F or C display, then release the button.				
	3.	Press the up and down arrow buttons O to toggle between the temperature scales.				
	4.	Press the standby button (b) to turn the screen on.				
Result	The tem	perature scale has now been changed.				

How to Change the Sound Settings

Before you begin	The warmer must be connected to electric power.						
Procedure	To change the sound settings, do the following.						
	Step	ep Action					
	1.	Press and while holding the temperature recall button , press the down arrow button . Once the sound volume setting displays, release both buttons.					
	2.	Press the up and down arrow buttons 🕢 妏 to adjust the volume. The volume range is 0 (mute) to 12 (maximum).					
	3.	Press the standby button 🕧 to turn the screen on.					
Result	The sou	ind setting has now been changed					

Result

The sound setting has now been changed.



How to Lock and Unlock the Controller

Before you begin	The warmer must be turned on (screen is on).					
Background	The controller can be locked to prevent changes being made to the temperature set-point.					
Locking the controller	To lock the controller, do the following.					
	Step Action					
	1.	Press and while holding the standby button (), press the up arrow button . The lock indicator in will illuminate.				
		NOTE: Make sure the set-point temperature did not change during the locking process.				
Unlocking the controller	To unloc	k the controller, do the following.				
	2.	Press and while holding the standby button , press the down arrow button . The lock indicator goes off.				
		NOTE: Make sure the set-point temperature did not change during the unlocking process.				
Result	The controller has now been locked or unlocked.					


How to Download Temperature Data from the Warmer (if equipped)

Before you begin	 The warmer must be turned on (screen is on). Do not remove the USB drive during the download process. You will need a SanDisk Cruzer Glide USB drive. 							
	NOTICE Only use a SanDisk Cruzer Glide USB drive larger than 4 GB. Using other USB drives may cause erratic operation or loss of data.							
Background	The warmer automatically downloads only the data that has not been previously downloaded. The warmer will store multiple years of data. The download time for one week of data should be under five seconds.							
	The warmer records various performance parameters during the following conditions:							
	automatically every five minutes while the power switch is ON.							
	each time the temperature set-point is changed.							
	each time the Top/TriMode cavity door is opened.							
	each time the fluid operation mode is changed.							
Procedure	To download temperature data from the warmer, do the following.							
	Step Action							
	1. Open the bottom cavity door.							

Continued on next page



Continued from previous page

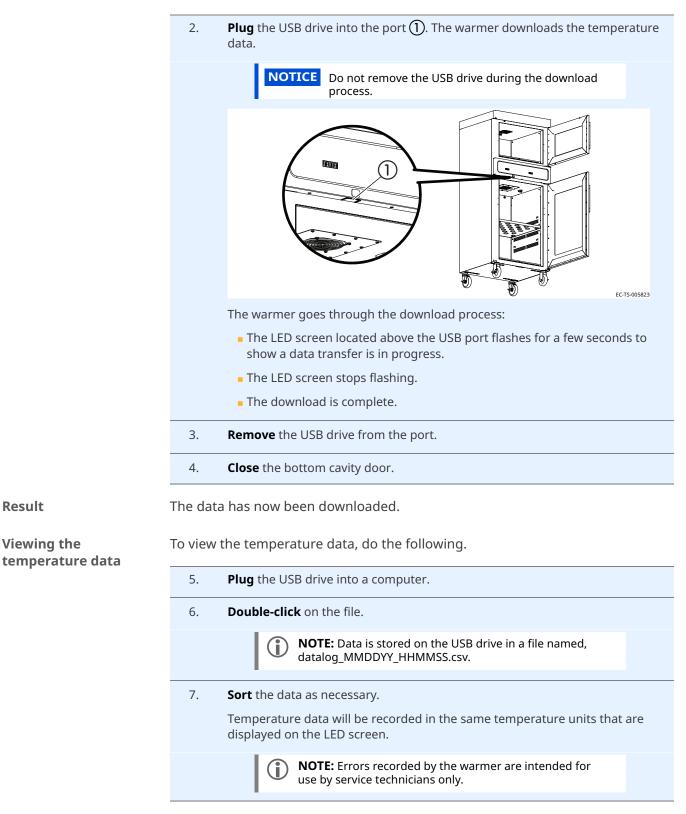

Continued on next page

Result

Viewing the



Continued from previous page

Viewing data in Microsoft® Excel The table below is an example of temperature data from the warmer.

Time Stamp	Chamber	Mode	Power	Door	Fan	Temperature	Set-point	Errors
12/19/2018 14:48	bottom		off	closed	off	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:48	top		off	closed	off	71.2	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:53	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:53	top	INJ	on	closed	on	70.4	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:58	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:58	top	INJ	on	closed	on	70.4	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 15:03	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 15:03	top	INJ	on	closed	on	70.3	104	Diagnostics.SensorSho rtErrors=0x0c

Result

The temperature data has now been viewed.



Maintenance Schedule

NOTE: Do not remove exterior panels from the warmer or attempt repairs. The warmer has no user-serviceable internal components. Only perform routine cleaning and maintenance procedures specifically described in this manual.

For daily maintenance, do the following.

Check:

- the air vents in the airflow insert panels are not obstructed (if applicable).
- all fan guards are clear and not obstructed (if applicable).
- the number of bottles/bags, or blankets as applicable do not exceed the maximum capacity per shelf or basket.

Clean:

any spills with a clean, lint free cloth. See topic *How to Clean the Warmer* for the list of approved cleaners.

Monthly

Daily

For monthly maintenance, do the following.

Check:

- the door gasket for tears and holes. Make sure that it is firmly attached to the door. Check the seal when the door is closed.
- the guards around the air temperature sensor are in place and fully secure to the warmer.
- the hardware securing the warmer(s) to the wall, if applicable.

Clean:

- the outside of the warmer with stainless steel cleaner. See topic *How to Clean the Warmer* for the list of approved cleaners.
- vacuum (if applicable):
 - \square fan openings
 - fan sail switch
 - vent openings



Yearly

For yearly maintenance, do the following.

Check:

- the set-point temperature compared to the actual temperature displayed.
 - Check the cavity air temperature with a thermocouple placed 1" (25mm) from the cavity sensor. Do not allow the thermocouple to touch any surface. Monitor the temperature for approximately one hour in an empty cavity.
- the condition of the plug and cord and replace if damaged.
- the controller screen for excessive wear. Make sure the controller screen displays and operates properly.

Clean:

the shelves and interior of the warmer. See topic *How to Clean the Warmer* for the list of approved cleaners.



How to Clean the Warmer

Before you begin



WARNING: Electric shock hazard. Disconnect the warmer from electric power before cleaning.



WARNING: Crushing hazard. An unstable warmer can lead to a crushing hazard. Observe

your facility's best practices for moving large equipment.



CAUTION: Burn hazard. Allow the warmer to cool before cleaning.

NOTICE Do not use:

- abrasive cleaning compounds.
- chloride based cleaners.
- commercial or household cleaners containing ammonia.

Approved cleaning agents

70% or greater isopropyl alcohol

PDI Sani-Cloth® HB

Monthly cleaning procedure

To clean the warmer monthly, do the following.

Step	Action
1.	Make sure that the warmer is disconnected from electric power and cool.
2.	Wipe the outside of the warmer with a stainless steel cleaner.
3.	Vacuum lint from fan openings, fan sail switch, and vent openings.
4.	Clean the interior of the warmer with a damp cloth or approved cleaning agent.
5.	Dry the interior and exterior of the warmer with a clean, lint-free cloth.
	Leave the door open until the interior of the warmer completely dries.
6.	Clean each side of the window pane with glass cleaner or distilled vinegar.




Yearly cleaning procedure

To clean the warmer yearly, do the following.

Step Action 1. Make sure that the warmer is disconnected from electric power and cool. 2. Remove the shelves ① from the interior of the warmer.



3. **Remove** the side panels (2) from the interior of the bottom cavity.

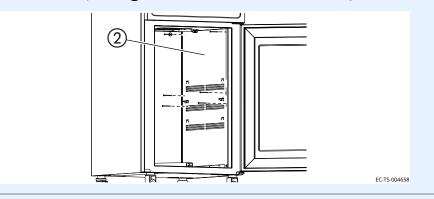

- 4. **Clean** the shelves and interior with a damp cloth or approved cleaning agent.
- 5. **Dry** the interior of the warmer with a clean, lint-free cloth.

Leave the door open until the interior of the warmer completely dries.

- 6. **Clean** each side of the window pane with glass cleaner or distilled vinegar.
- 7. **Re-install** the side panels and shelves.



WARNING: Personal injury hazard.

The blankets or fluids may overheat if the side panels and shelves are not in place. Make sure to re-install the side panels and shelves before operating the warmer.

Result

The warmer is now clean.



Warmer Disposal / Decommissioning

This product and its accessories must be disposed of according to local laws and regulations. Do not dispose of this product as unsorted municipal waste.



This page is intentionally left blank.



What to do if a Power Interruption Occurs

Background	You may need to reset the warmer in the case that a power interruption occurs. When the power is restored, the ON/OFF status indicator LED decimal flashes LED The controller stores all settings and continues operation using these settings in the event of a power interruption.	
Procedure	To contir	nue operation of the warmer, do the following.
	Step	Action
	1.	Press the standby button (). The status indicator LED decimal goes out. If an error code displays, refer to topic <i>Error Codes</i> .
	2.	Make sure the fluids or blankets are at the correct temperature before using.
		NOTE: Inspect the fluids for their expiration date before use.
Result	Resume operation of the warmer.	



What to do if the Alarm Indicator Light Flashes

Background	The alarm indicator light of flashes and the controller sounds an alarm when the warmer malfunctions.		
Procedure	If the alarm indicator light flashes, do the following.		
	Step	Action	
	1.	Press the standby button () to acknowledge and mute the alarm.	
	The alarm indicator light goes off. Refer to the troubleshooting section for the error code.		
	2.	Make sure the fluids or blankets are the correct temperature before using.	
Result	Resume operation of the warmer.		



Error Codes

Background

This section is provided for the assistance of qualified and trained service technicians only and is not intended for use by untrained or unauthorized service personnel. Failure to observe this precaution may void the warranty.

NOTE: If the warmer is not operating properly, check the following before calling an authorized service agent:

Verify that the power to the warmer is on.

If applicable, ensure the female end of plug is securely seated in the warmer and that the male end of plug is in an appropriate, functioning outlet.

If applicable, examine the fuses. Replace the fuses, refer to topic *How to Replace a Fuse*.

If applicable, examine the high limit manual reset button. If the high limit manual reset button is tripped, reset the warmer, refer to topic *How to Manually Reset the Warmer*.

NOTE: All non-critical codes can be cleared using the standby button. Critical errors (marked with a *) can only be cleared by setting the power switch at the rear of the warmer to the off (O) position and allowing the warmer to cool.

NOTICE Do not attempt to repair or service the warmer beyond this point. Contact the manufacturer for the nearest authorized service agent. Repairs made by any other service agent without prior authorization by the manufacturer will void the warranty.

Code	Refers to	Action required
Display flashes set point	Cavity temperature higher than set point	Cavity temperature is higher than the set point temperature. Allow cavity to cool to set point temperature.
E-10 ES10 ES20 ES30 ES40 ES50 ES60 ES70	Cavity sensor Sensor 1 Sensor 2 Sensor 3 Sensor 4 Sensor 5 Sensor 6 Sensor 7	 Sensor is shorted. Software disengages heating pads. Acknowledge error by pressing the standby button. If error persists, a qualified service technician should test the sensor. Test the sensor. Detach the sensor from the warmer. Use an Ohm meter to measure the resistance of the sensor. Check the sensor at 25°C (77°F). If the reading is 10 KOhm ±1.5 KOhm, replace the display. If the reading is ±2 KOhm, replace the sensor. Check the wires for integrity. Inspect the connections at the control and terminal block to ensure proper and secure connections. If necessary, re-secure the faulty connections. Contact service if error persists.

TROUBLESHOOTING

Code	Refers to	Action required
E-11 ES11 ES21	Cavity sensor Pad sensor 1 Pad sensor 2	Sensor is open. Software disengages heating pads. Acknowledge error by pressing the standby button. If error persists, a qualified service technician should test the sensor.
ES31 ES41 ES51 ES61 ES71	Pad sensor 3 Pad sensor 4 Pad sensor 5 Pad sensor 6 Pad sensor 7	Test the sensor. Disconnect the sensor from the warmer control PCB and use an Ohm meter to measure the resistance of the sensor. Sensor PR-37140 should measure 10 KOhm ± 1.5 KOhm at 25°C (77°F). Sensor SN-33541 should measure 100 Ohm ± 10 Ohm at 25°C (77°F). Replace sensor if measurement is outside of tolerance.
		Check the wires for integrity. Check for proper and secure connections at the control and terminal block. If necessary, re-secure the faulty connections.
		Contact service if error persists.
E-30	Under temperature	The cavity temperature is lower than the set temperature for 90 minutes or longer.
		■ Make sure the door is closed.
		If the cavity is overloaded, redistribute the inventory. Do not exceed the height of the insert.
		Test the sensor. Disconnect the sensor from the warmer control PCB and use an Ohm meter to measure the resistance of the sensor. Sensor PR-37140 should measure 10 KOhm ± 1.5 KOhm at 25°C (77°F). Sensor SN-33541 should measure 100 Ohm ± 10 Ohm at 25°C (77°F). Replace sensor if measurement is outside of tolerance.
*E-31	Cavity sensor	The sensor reading is above the temperature set-point. Blanket warming mode triggers at 15° over set-point. Fluid warming mode triggers at 5° over set-point.
		The difference between the room temperature and the fluid set-point temperature must be greater than 11°C (20°F).
		Contact service if error persists.
P131	Pad sensor 1	Heater pad over-temp error.
P231 P331	Pad sensor 2 Pad sensor 3	■ Software disengages heating pads.
P431	Pad sensor 4	Acknowledge error by pressing the standby button.
P531 P631	Pad sensor 5 Pad sensor 6	Allow the warmer to cool.
P731	Pad sensor 7	Contact service if the error persists.
*E-33	Cavity sensor	 Sensor reading is above maximum allowable temperature set-point and over temp value. Blanket warmers trigger at 82°C (180°F). Fluid warmers trigger at 71°C (160°F).
		Contact service if error persists.
*E-50	Analog to Digital Converter Error	Remove inventory, discard if necessary, and allow the warmer to cool down.
		If error persists after cool down and reset, the control assembly should be replaced by a qualified service technician. Contact service.



Code	Refers to	Action required	
E-60	Real Time Clock	Real time clock rechargeable battery backup has discharged.	
	Checksum Error (Blanket warmers only)	Plug the warmer into the outlet for 30 minutes.	
*E-61	Real Time Clock (Blanket warmers only)	Real time clock not responding. Contact service if error persists.	
E-62	Real Time Clock (Blanket warmers only)	Timer overlay is present, but no real time clock is detected. Contact service.	
*E-70	Pad Count Error	More heater pads detected than set for. Hold the standby button for 12 seconds until display shows "PAd#" (# = number of pads selected [3-7]). Press up or down arrow to adjust to correct number of pads. [Blanket warmers: 2.5ft ³ = 2 pads, 3.5ft ³ , 4.0ft ³ , 7.5ft ³ = 4 pads.	
		Fluid warmers: 2.5ft ³ & 4.0ft ³ = 3 pads].	
*E-71	Personality Error	Contact service.	
E-80	EEPROM Error	EEPROM not responding. Contact service if error persists.	
*E-81	Calibration not locked	Contact service.	
*E-83	EEPROM Error	Contact service for help resetting the control.	
E-87	EEPROM Error	Stored offsets corrupted. Offsets reset to 0.	
		Control may need to be re-calibrated.	
		Possible bad EEPROM.	
		Contact service if error persists.	
E-90	Button stuck	A button has been held down for >60 seconds. Adjust control. Error will reset when the problem has been resolved.	
E-95	Factory Test pin short detected.	Ensure that debris is not causing a short between the test pins. If the pins are clean, replace the control.	
*E-98	Temperature Delta Error	Temperature of the cavity sensors 1 and 2 differ by more than 3.3°C (6°F). Remove product and allow the warmer to cool down.	
		Verify that the product sensor is clean and operating correctly.	
		Set the power switch to the OFF position to clear the error code.	
		If error persists, the cavity sensor should be replaced by a qualified service technician. Contact service.	
E-99	Hardware Over Temp	Inspect the connections and condition of high limit bimetal thermostat.	
		If applicable, make sure that the compartment fan motor is operating.	
		Contact service if error persists.	

TROUBLESHOOTING

Code	Refers to	Action required
*EFAn	Fan or Fan Sensor Failure	 Check to make sure the fan sensor wires did not disconnect from the fan or the control board. If the error persists after checking the wires, replace the fan.
Flashing yellow iv fluid temperature status indicator	Flashing yellow LED is illuminated within iv storage shelf	 Set the power switch to the OFF position. Remove any iv bags from inside the warmer. Set the power switch to the ON position. Visually check to see if flashing yellow LED remains illuminated. Contact service if error persists.

*All non-critical codes can be cleared using the standby button. Critical errors (marked with a *) can only be cleared by setting the power switch at the rear of the warmer to the OFF position and allowing the warmer to cool.



Guidance and Manufacturer's Declaration

The warmer requires special precautions regarding EMC (Electromagnetic Compatibility) and needs to be installed and put into service according to the EMC information provided in the accompanying documents.

Portable and mobile RF communications equipment can affect medical electrical equipment.

A risk of increased emissions or decreased immunity may result if the power cord is altered or a manufacturer supplied power cord is not used.

The warmer should not be used adjacent to or stacked with other equipment.

The essential performance of the warmer is to not exceed an internal temperature of 160°F (71°C) for blanket warmers or 150°F (66°C) for fluid warmers.

Electromagnetic emissions

The warmer is intended for use in the electromagnetic environment specified below.

Emission test	Compliance	Electromagnetic environment - guidance	
RF emissions; CISPR 11	Group 1	The warmer uses RF energy only for internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.	
RF emissions; CISPR 11	Class B	The warmer is suitable for use in all establishments, including domestic	
Harmonic emissions; IEC 61000-3-2	N/A	establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.	
Voltage fluctuations/Flicker emissions; IEC 61000-3-3	N/A		



ELECTRICAL INFORMATION

Continued from previous page

ElectromagneticThe warmer is intended for use in the electromagnetic environment specifiedimmunitybelow.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electromagnetic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines; ±1 kV for input/output lines	+2 kV for power supply lines	Main power quality should be that of a typical commercial or hospital environment. The warmer does not have any input/output lines.
Surge IEC 61000-4-5	±1 kV differential mode; ±2 kV common mode	±1 kV differential mode; ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT (>95% dip in UT) for 0.5 cycle 40% UT (60% dip in UT) for 1 cycle 70% UT (30% dip in UT) for 25 cycles <5% UT (>95% dip in UT) for 250 sec	<5% UT (>95% dip in UT) for 0.5 cycle 40% UT (60% dip in UT) for 1 cycle 70% UT (30% dip in UT) for 25 cycles <5% UT (>95% dip in UT) for 250 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the warmer requires continued operation during power mains interruptions, it is recommended that the warmer be powered from an uninterrupted power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: UT is the a.c. mains voltage prior to application of the test level.



ElectromagneticThe warmer is intended for use in the electromagnetic environment specifiedemissionsbelow.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 V/m 150 kHz to 80 MHz	3 V/m	Portable and mobile RF communications equipment should be used no closer to any
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.7 GHz		part of the warmer, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
			Recommended separation distance
			d = [3.5/3] √P
			d = [3.5/3] √P 80 MHz to 800 MHz
			d = [7/3] √P 800 MHz to 2.5 GHz
			where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).
			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ¹ , should be less than the compliance level in each frequency range ² .
			Interference may occur in the vicinity of equipment marked with the following symbol:

NOTE: At 80 MHz and 800 MHz, the higher frequency range applies. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

1. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the warmer is used exceeds the applicable RF compliance level above, the warmer should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the warmer.

2. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [VI] V/m.



Electromagnetic immunity distance

The warmer is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The user of the warmer can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the warmer as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = \left[\frac{3.5}{3}\right] \sqrt{P}$	$d = \left[\frac{3.5}{3}\right] \sqrt{P}$	$\mathbf{d} = \begin{bmatrix} \frac{7}{3} \end{bmatrix} \sqrt{\mathbf{P}}$
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.738
1	1.167	1.167	2.333
10	3.689	3.389	7.379
100	11.667	11.667	23.333

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.



Warranty

Introduction	Pedigo Products, Inc. warrants to the original purchaser that any original part that is found to be defective in material or workmanship will, at our option, subject to provisions hereinafter stated, be replaced with a new or rebuilt part.		
Warranty Period	The original parts warranty period is as follows:		
	The labor warranty remains in effect one (1) year from installation or fifteen (15) months from the shipping date, whichever occurs first.		
	The original parts warranty for the cavity fan motor remains in effect one (1) year from installation of appliance or fifteen (15) months from the shipping date, whichever occurs first. The original parts warranty on all other parts remains in effect three (3) years from installation of appliance or thirty-nine (39) months from the shipping date, whichever occurs first.		
Exclusions	This warranty does not apply to:		
	Calibration		
	 Equipment damage caused by accident, shipping, improper installation or alteration. 		
	 Equipment used under conditions of abuse, misuse, carelessness or abnormal conditions including equipment subjected to harsh or inappropriate chemicals including but not limited to compounds containing chloride or quaternary salts, poor water quality, or equipment with missing or altered serial numbers. 		
	 Any losses or damage resulting from malfunction, including loss of contents or consequential or incidental damages of any kind. 		
	 Equipment modified in any manner from original model, substitution of parts other than factory authorized parts, removal of any parts including legs, or addition of any parts. 		
	Collateral or incidental damage as a direct result of servicing equipment built into a wall structure is not covered under warranty. It is the responsibility of the owner to bear all expense related to structural repairs including, but not limited to, external electrical connections and wiring, and the removal or replacement of caulk, grout, tile, or wall covering of any kind. A service access panel for built-in equipment installations is strongly recommended.		
Conclusion	This warranty is exclusive and is in lieu of all other warranties, expressed or implied, including the implied warranties of merchantability and fitness for purpose. In no event shall the company be liable for loss of use, loss of revenue, or loss of contents or revenue, or for indirect or consequential damages. This warranty is in lieu of all other warranties expressed or implied and Pedigo Products, Inc. neither assumes or authorizes any persons to assume for it any other obligation or liability in connection with Pedigo Products, Inc. equipment.		





Pedigo Products, Inc. 4000 S.E. Columbia Way Vancouver, WA 98661 Phone: 800.822.3501 | Fax: 360.696.1700

www.pedigo-usa.com

Printed in the U.S.A.

Specifications are subject to change without notice.

Made in the U.S.A